

# Název:Metalothionein a jeho vztah k<br/>rakoviněŠkolitel:Vojtěch AdamDatum:31.1.2013

Reg.č.projektu: CZ.1.07/2.4.00/31.0023

Název projektu: Partnerská síť centra excelentního bionanotechnologického výzkumu



#### Content

I. What are metallothioneins?

II. The biologically important roles of metallothioneins


#### Content

I. What are metallothioneins?

II. The biologically important roles of metallothioneins

#### **Metallothioneins = proteins**

- Intracellular, low molecular and cysteine-rich proteins with molecular weight from 6 to 10 kDa
- MTs consist of two binding domains  $-\alpha$  and  $\beta$ .
- N-terminal part of protein  $-\beta$ -domain; three binding places for divalent ions.
- C-terminal part of protein  $\alpha$ -domain; four binding places for divalent ions.
- The most repeated structure motif: cysteine(C)–serine(S)–cysteine(C).

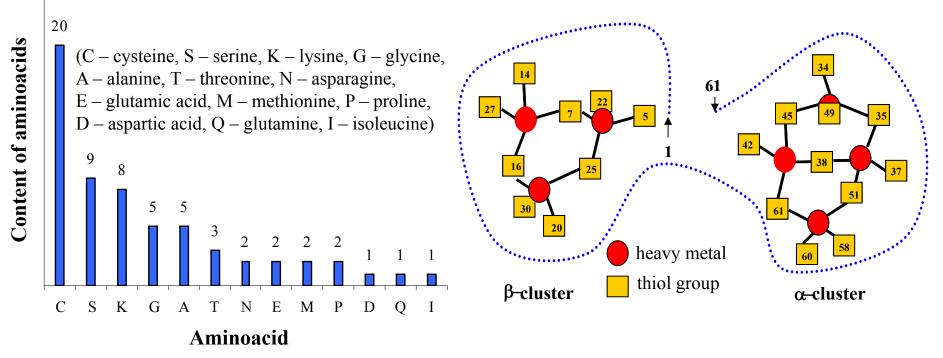


Adopted from J. Petrlova, et al. Attomole voltammetric determination of metallothionein, Electrochim. Acta 51 (2006) 5112-5119.

#### Content

I. What are metallothioneins?


## II. The biologically important roles of metallothioneins



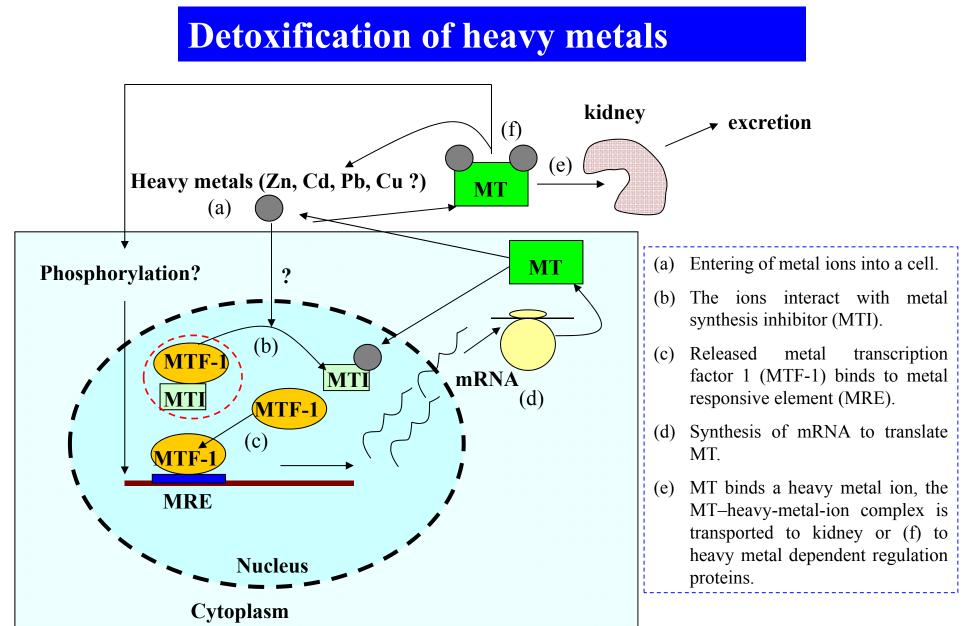


Adopted from Eckschlager, et al. (2009). Metallothioneins and Cancer. Curr. Protein Pept. Sci., 10, 360-375.

#### **Metallothioneins and NF-ĸB?**

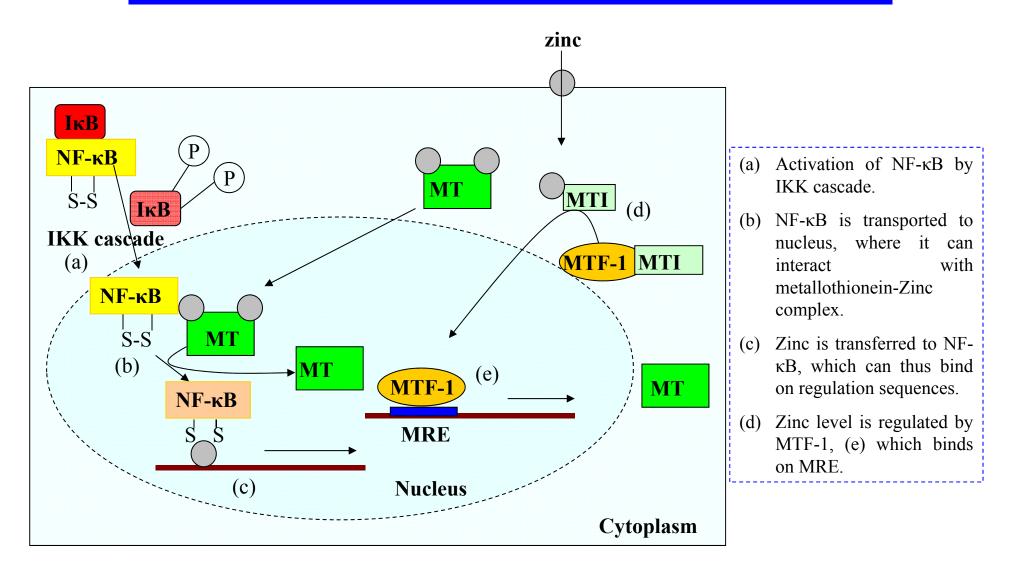



Adopted from Eckschlager, et al. (2009). Metallothioneins and Cancer. Curr. Protein Pept. Sci., 10, 360-375.

#### Metallothioneins as scavengers of reactive oxygen species

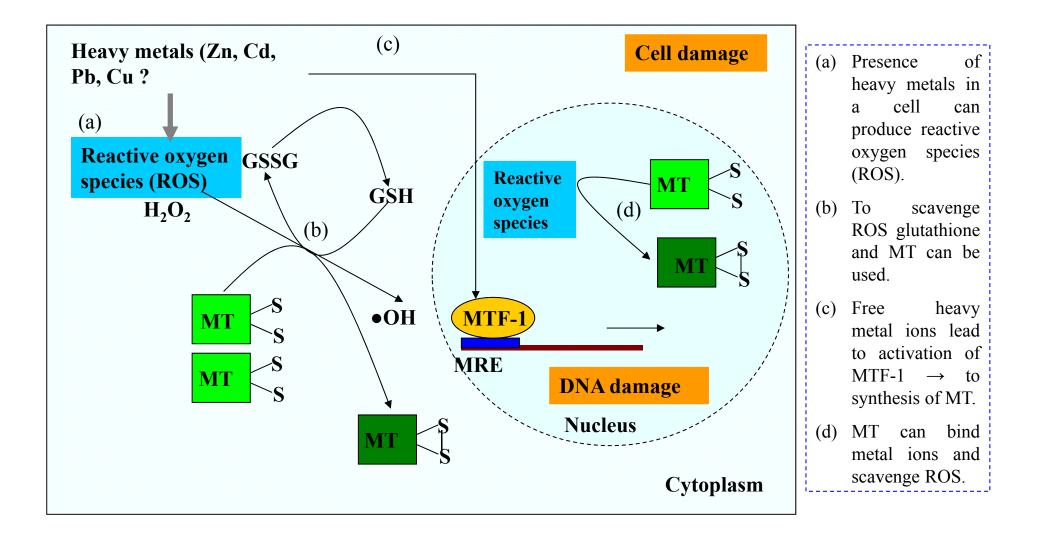


Adopted from Eckschlager, et al. (2009). Metallothioneins and Cancer. Curr. Protein Pept. Sci., 10, 360-375.

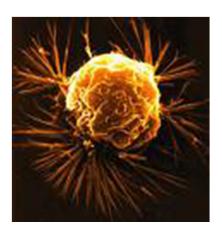


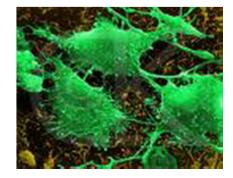



#### **Can metallothioneins play a key role in cancerogenesis?**

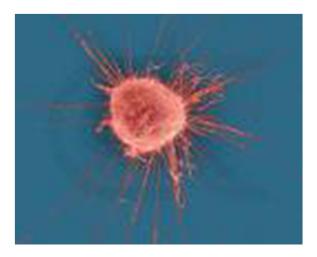


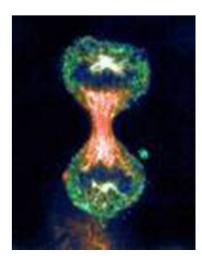

#### Content

I. What are metallothioneins?


II. The biologically important roles of metallothioneins

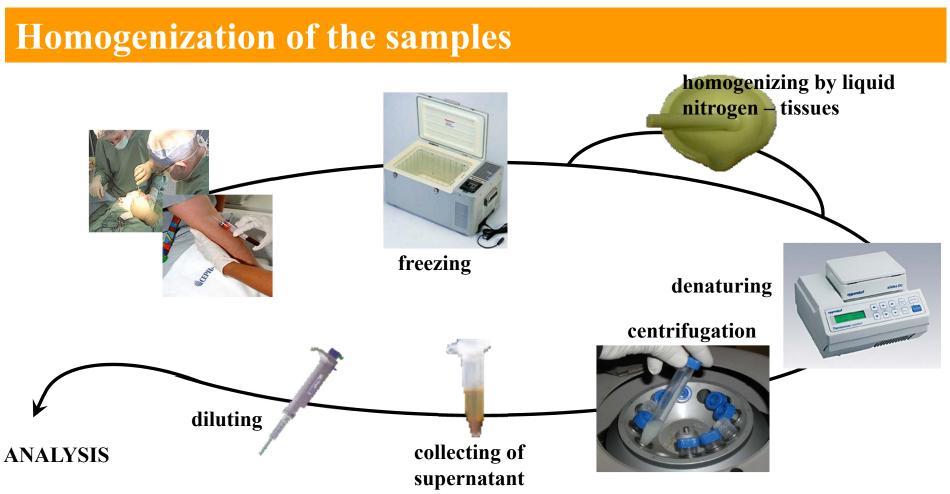


- Briefly, the sample was kept at 99 °C for 15 min. with occasional stirring, and then cooled to 4 °C. The denatured homogenates were centrifuged at 4 °C, 15 000 g for 30 min. Heat treatment and solvent precipitation effectively denature and remove high molecular weight proteins out from samples; metallothionein belongs to thermo stable proteins.
- The prepared samples are analysed by Adsorptive Transfer Stripping Technique coupled with Differential Pulse Voltammetry Brdicka Reaction.

#### **Brdicka reaction**

- Brdicka reaction the hydrogen evolution from supporting electrolyte containing 1 mM  $Co(NH_3)_6Cl_3$  and 1 M ammonia buffer  $(NH_3(aq) + NH_4Cl, pH = 9.6)$  in the presence of peptides and/or proteins containing thiol group.
- limit of quantification: **50 pM** (1 fmol MT in 5  $\mu$ l).

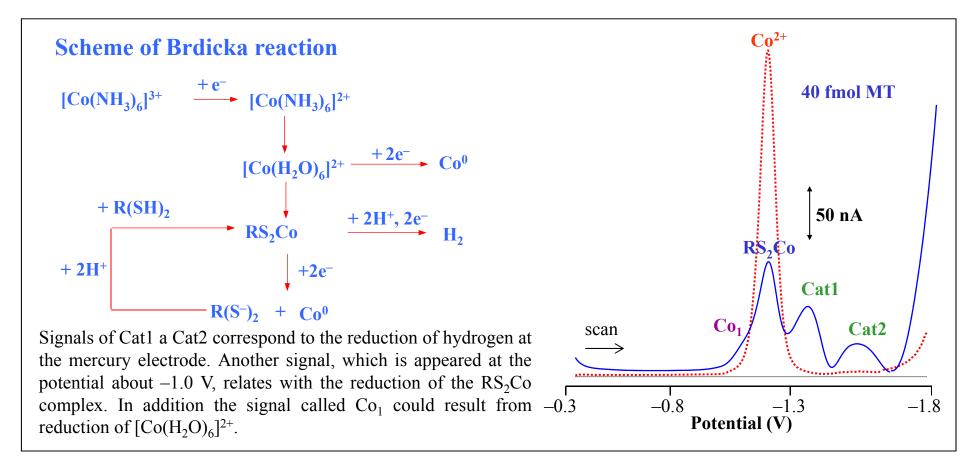

Adapted from J. Petrlova, et al. Attomole voltammetric determination of metallothionein, Electrochim. Acta 51 (2006) 5112-5119.

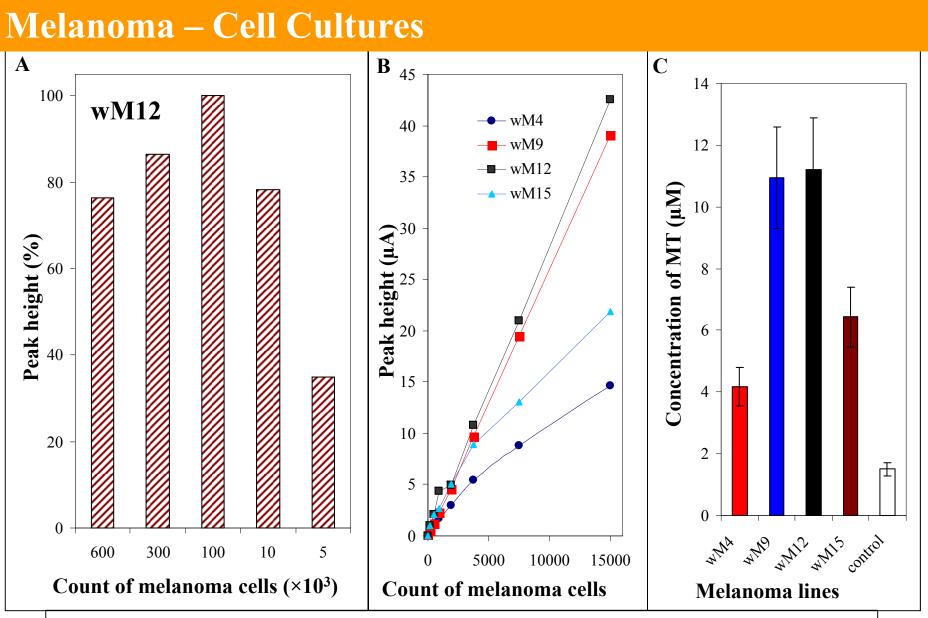

*MT* level determined at control cells without symptoms of tumour transformation  $-1.5 \mu M$ .

Adopted from S. Krizkova, et al. Utilizing of adsorptive transfer stripping technique Brdicka reaction..., Sensors 8 (2008) 3106-3122...


#### Melanoma – Tissue

- Using directed selection an original cancer model was established in the Institute of Animal Physiology and Genetics in Liběchov a strain of miniature pigs that was designated with acronym MeLiM (Melanoma-bearing Libechov Minipig).
- Melanoma in this strain is heritable.
- Multiple skin nodular tumours (i.e. the most aggressive form of melanoma) appear on various parts of body in about a half of piglets.
- Their histological, immunohistochemical and biochemical characterization and a broad melanoma cell dissemination document similarities with human melanoma and malignant behaviour of this porcine cancer.




Adopted from S. Krizkova, et al. Utilizing of adsorptive transfer stripping technique Brdicka reaction..., Sensors 8 (2008) 3106-3122...

#### Melanoma – Tissue

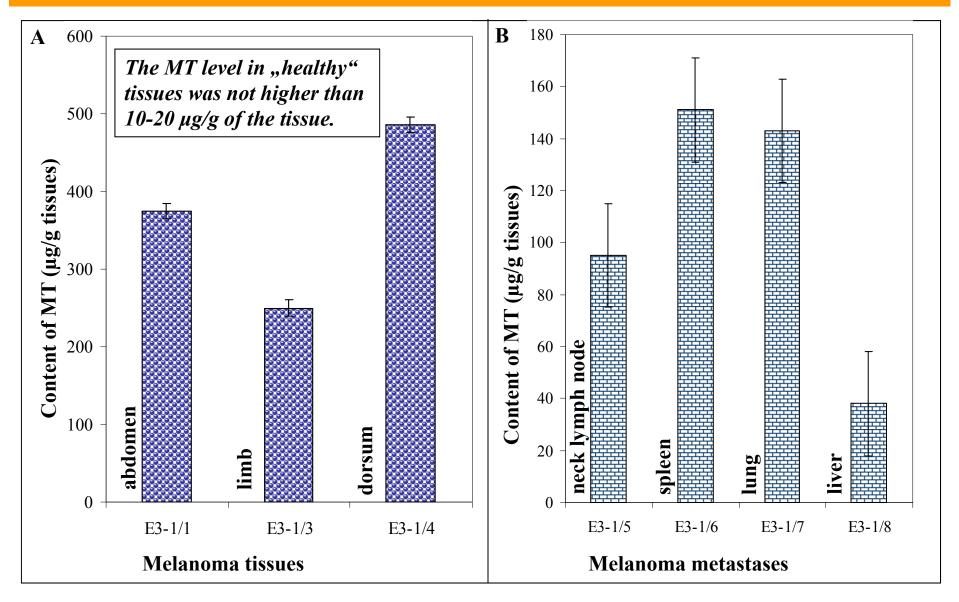

Adopted from S. Krizkova, et al. Utilizing of adsorptive transfer stripping technique Brdicka reaction..., Sensors 8 (2008) 3106-3122...

#### **Blood Serum Samples from Patients with Melanoma**

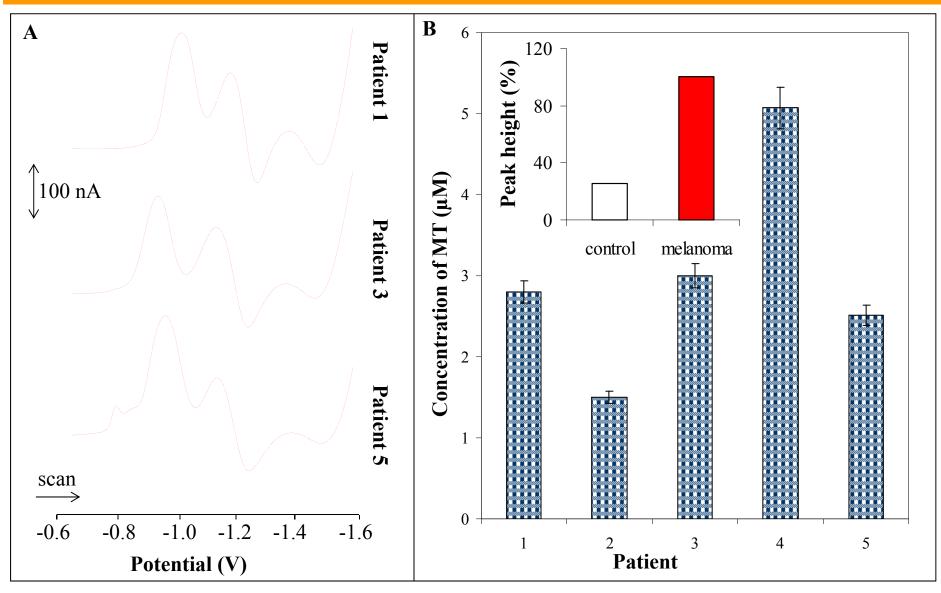

Adopted from S. Krizkova, et al. Utilizing of adsorptive transfer stripping technique Brdicka reaction..., Sensors 8 (2008) 3106-3122...



### Thank you very much

Reg.č.projektu: CZ.1.07/2.4.00/31.0023

Název projektu: Partnerská síť centra excelentního bionanotechnologického výzkumu

